

MID SKIN SOLUTIONS
LABORATORIES

Pluryal
mesoline system

Pluryal mesoline system

Your Needle-Free Mesotherapy Solution

#bemesoline


your most beautiful age is now
P

1. Experience Needle-Free Mesotherapy with the Mesoline System

Effective	Glowing Results	Painless
No Downtime	Social Friendly	Needle-free

2. About the Device

MD Skin Solutions, creators of the trusted Mesoline sterile serums, have teamed up with Italy's Panestetic to create a **groundbreaking needle-free mesotherapy machine**. This **fast, efficient, and painless** device offers beauty professionals a **complete solution for radiant, glowing skin**.

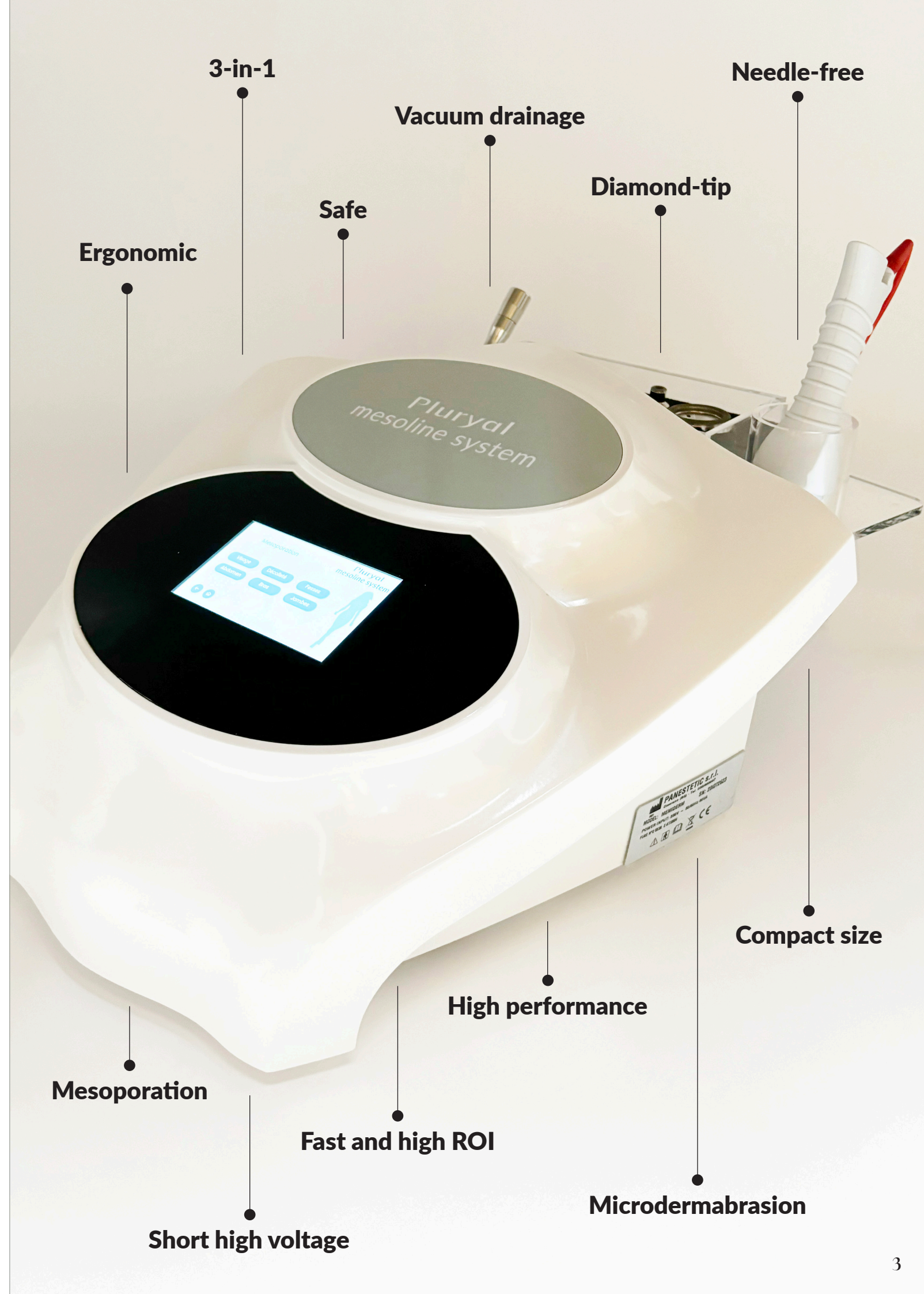



3. The Mesoline System: A Three Step Treatment

I. Vacuum Drainage

A mechanical massage that **applies negative pressure** to the skin, **enhancing microcirculation** and **toxin removal**.

Benefits:

- **Reduces** puffiness, bags, and dark circles
- **Brightens** skin tone
- **Boosts** regeneration and circulation

II. Microdermabrasion

Diamond-tip exfoliation for treating uneven skin tone, scars, and pigmentation. **Superior to other types of dermabrasion** and ideal for all skin types.

Benefits:

- **Removes** dead skin cells and impurities
- **Reduces** pigmentation and wrinkles
- **Soothes and brightens** the skin



III. Mesoporation

The stratum corneum is the outer layer of the skin which consists of keratin-filled cells surrounded by lipid bilayers, which restrict the movement of substances across the skin.

The mesoporation is a technique that enhances the skin's permeability to bypass the stratum corneum barrier effectively and then deliver active ingredients into the skin without needles for enhanced results.

4. How mesoporation works:

- Uses **short high-voltage pulses**
- **Increases the delivery rates** by up to 10,000 times
- **Rearranges the lipid bilayer structure of the stratum corneum** and applies an electrophoretic force
- After applying a pulse, **the membrane charges, becomes unstable, and electroporation occurs**, reducing membrane resistance, which is reversible
- **Increases the transport of both small ions and large molecules** across membranes

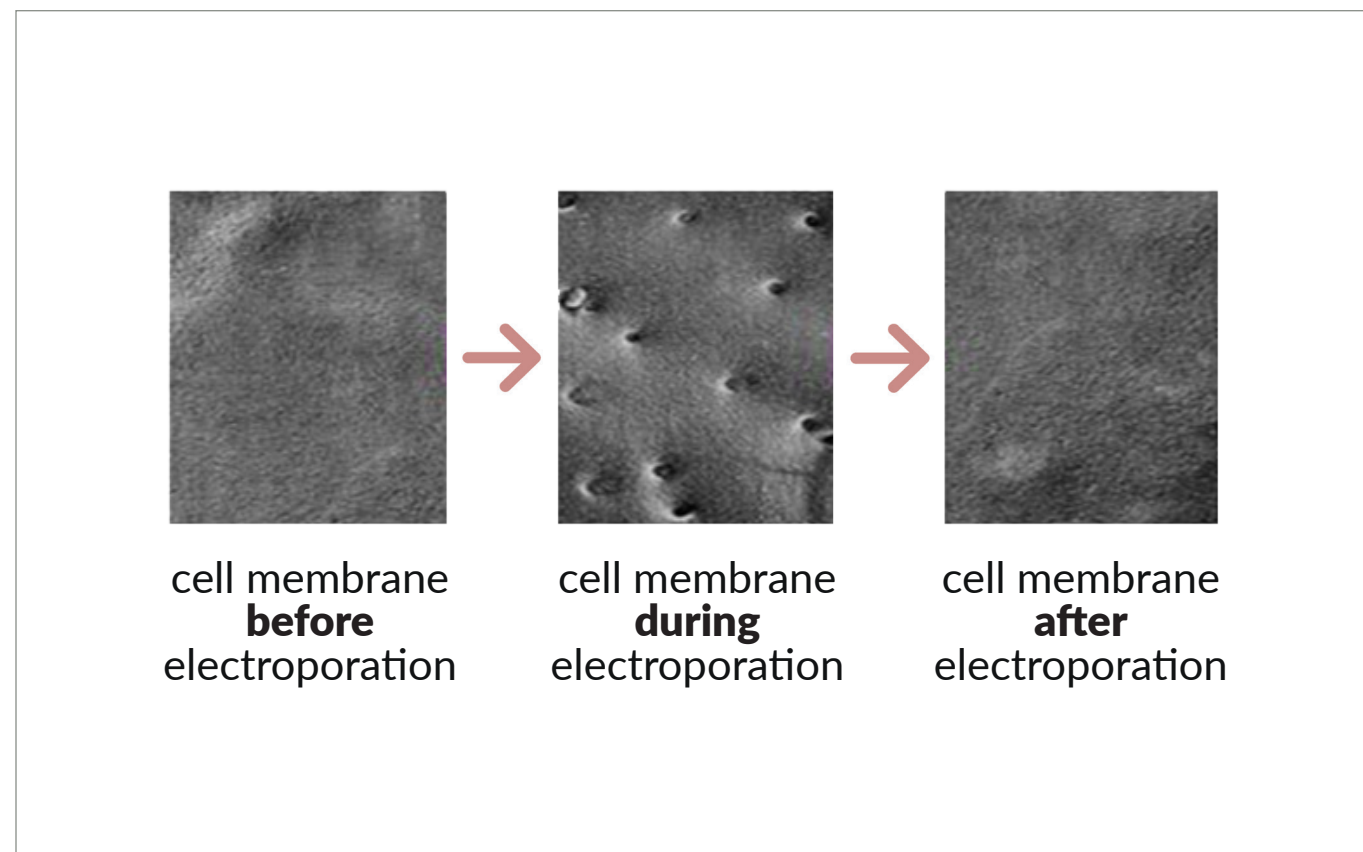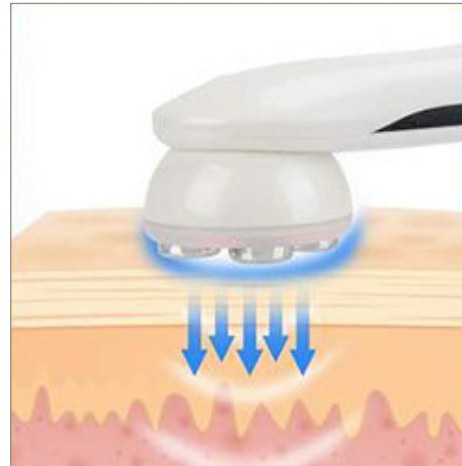

Face:


Delivers Mesoline treatments to reach the desired results:

- **Tightening, revitalizing, renewal, brightening, antioxidant**, etc.
- Improves the skin's texture by **reducing redness, melanin spots, and imperfections**
- **Improves the levels of hydration**
- **Reduces wrinkles**

Body:

Delivers Mesoline treatments to reach the desired results:

- **Attenuates cellulite**
- **Lifts and tones** the body
- **Enhances** microcirculation

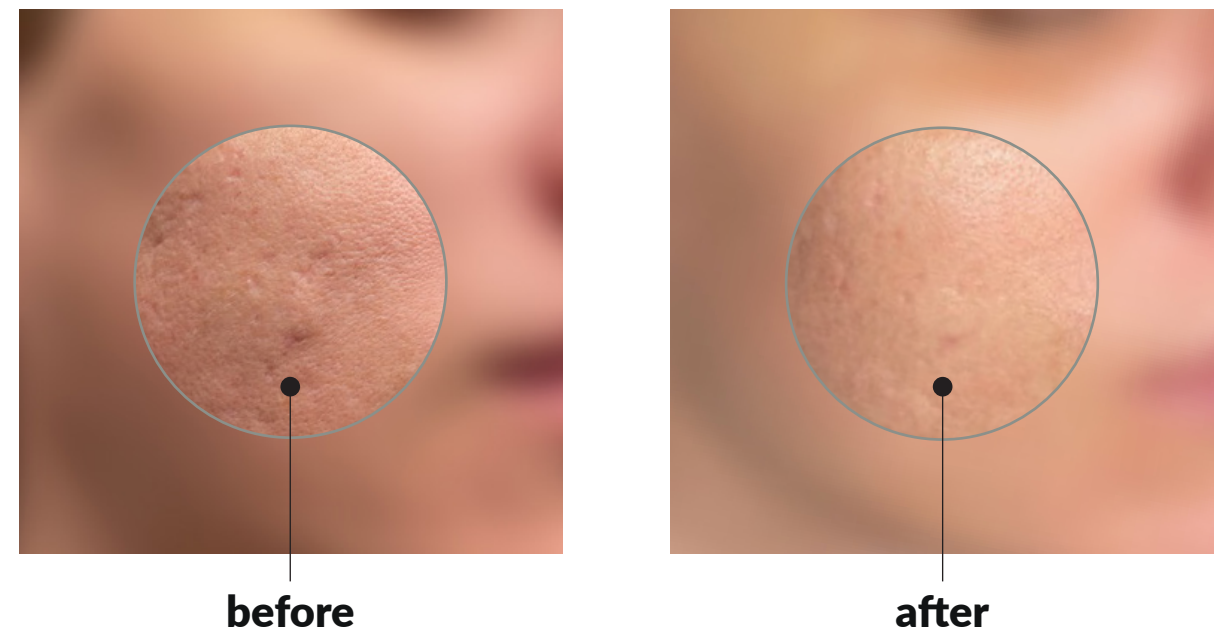

Courtesy of Dr. Marina Disteldorf, Luxembourg.

After 3 sessions:
Vacuum + microdermabrasion + mesoporation with Mesoline Refresh +
mask HA and 10 min of LED



Face:

Protocol	Products used	Dosage	No. of sessions	Frequency
Regenerating	Refresh + Exoyouth	3mL of each product	5-6 sessions	Every 2-3 weeks
Firming	Refresh + Tight			
Detoxifying	Antiox + Clear			
Clarifying	Refresh + Radiance			
Brightening	Refresh + Shine			

5. Protocols and Association with Mesoline

The Mesoline System is designed to be associated with the **Mesoline sterile serums.**



Body:

Protocol	Products used	Dosage	No. of sessions	Frequency
Firming / Cellulite	Bodyfirm	5-15mL	6-8 sessions	Every 2-3 weeks
Scars or stretch marks	Tight & Exoyouth			
Draining / Legs	Bodyfirm			
Hands	Refresh + Radiance (hyperpigmentation) or Tight (wrinkles or firmness)			
Neck and décolleté	Refresh + Radiance (hyperpigmentation) or Tight (wrinkles or firmness)			

Important: Full of treatment possibilities, Mesoline Serums can be used alone or combined → depending on the customer's skin needs.




• Vacuum and microdermabrasion handpiece



• Mesoporation handpiece with spheres



• 9 diamond tips of different strengths and sizes

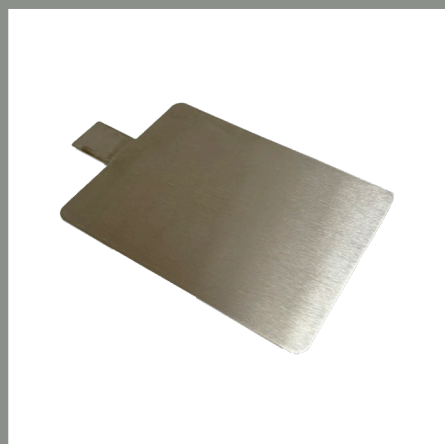

• Steel plate

6. Technical Data

Description :

Supply

Absorption

Fuses

Vacuum

Max power EP

Visualization parameters

Setting of the parameters

Commands

Safety

Controls

Alarms

Classification according to the Standard CEI 62-39

Operating temperature

Operating humidity

Weight

Languages

Characteristics :

230V ± 10% - 50/60Hz

50/60 VA (max)

2 x 2 A - type T - 5x20 (IEC 60127)

Adjustable from 0 to -800 millibar

30 Vpp ohm

Colour graphic display 4,3"

Icon menu and text

TOUCH SCREEN

Isolated power transformer a 4.000 V
Internal fuses

Microprocessor control of all technical and treatment parameters.
Feedback circuits on all outputs.

Fuses on hdmi board

Equipment class I
Appliance with parts to apply type BF
Device with casing not protected against the penetration of liquids (degree fi protection IP20)

From 10 to 40°C

From 30 to 75%

5 kg

French, English, German, Portuguese, Italian, Spanish, Czech, Finnish, Polish

Pluryal
mesoline system

How to order:

**Contact your local distributor for more information
or visit our website and social media pages:**

www.pluryal.com

[@pluryal](https://www.instagram.com/pluryal)

Contact:

MD SKIN SOLUTIONS
LABORATORIES

9AB boulevard du Prince Henri
L-1724 Luxembourg

New address from 2025 March, 1:

48, rue de Bragance
L-1255 Luxembourg

+352 26 20 04 06

info@mdskin-solutions.com